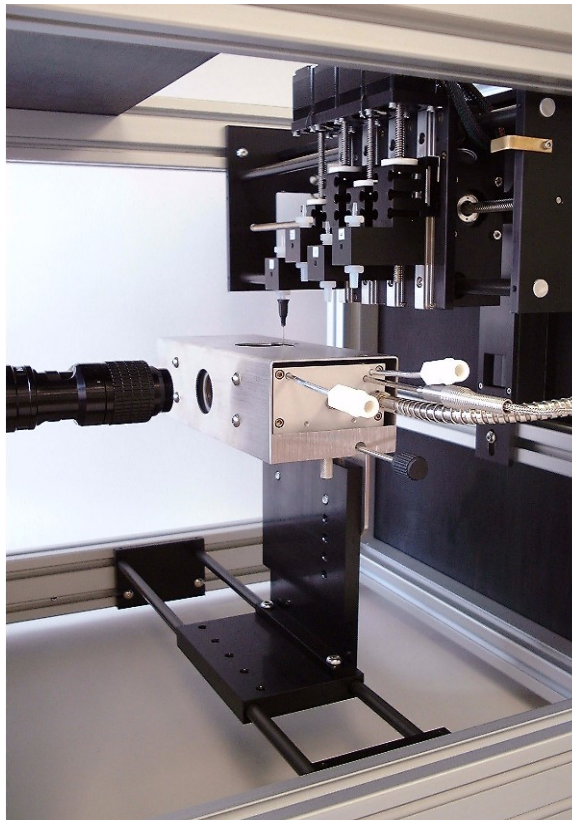

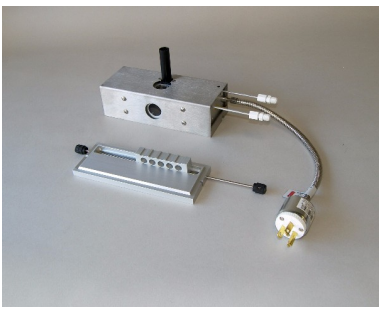
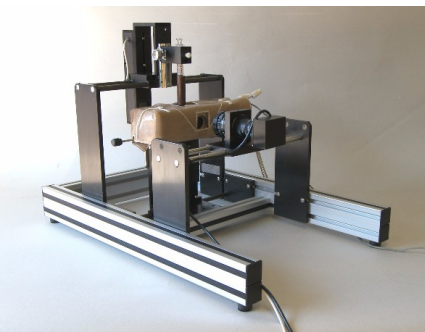


300°C Environmental Chamber

- Electrically heated chamber
- Contact angle measurements at elevated temperature
- Two gas ports are provided for purge gas
- Internal loop through body for user supplied circulating bath
- Two 20mm diameter windows provide light passage
- Top window allows sample inspection
- Specimen stage with manual linear slide holds five 10 x 15mm samples
- Needle port provided above sample
- Removable Kevlar cover to improve temperature stability and protect user from hot surface
- Includes stand-alone temperature controller
- Mounting hardware for FTA1000, FTA200 and FTA188



Controller	Chamber & Specimen Slide	Chamber Mounted on FTA1000B
		

For more information contact sales@firsttenangstroms.com